



Storage

1. Please cut off the power when not in use;
2. Do not dip water or other liquid into the product;
3. Please clean the product with a soft cloth dipped with water or neutral detergent.
4. Do not clean the product with water, remove this sign and just put.
5. Store the product in a dry, cool place, do not expose the product directly to sunlight, high temperature, fog or damp environment, remove this sign and just put.

Troubleshoot

Symptom	Cause/ Treatment
Unable to turn on the power	Please confirm that the massage button has been pressed.
Suddenly stops working during massage	1. Accidentally turned the switch off. 2. Product has 15 minute automatic timing device, time limit has been reached.

Technical Parameters

Product name: Dr.Physio Shiatsu Massager
Rated voltage: DC 12V Rated power: 24W

For any query Contact: customer@nureca.com



Shiatsu Massager



Instruction Manual



Please read the instructions carefully given in the manual.

ATTENTION


1. It is advised to consult a doctor for the pregnant women and patient with serious disease like malignancy.
2. Please do not touch the internal parts in any case.
3. Do not clean the switch with water.
4. Unplug the power cord when not in use.
5. Do not use a damaged power cord or plug.
6. Do not spray volatile oil, paint or pesticides.
7. Do not use in humid places such as bathrooms.
8. Please avoid direct sunlight for long periods of time to avoid damage.
9. Please do not open the shell or replace the internal parts, avoid short circuit and electric shock.
10. Do not place near fire.
11. Please keep this product away from children.

Operation instructions

1. Insert the plug into the charging jack.
2. Pressing "O" button, will open massage and warm function automatically.
3. Pressing "Q" button, will make the massage head turn clockwise or counter clockwise.
4. The product have three- speeds intensity to adjust, pressing "H" can adjust massager speed from low, medium and high.
5. Pressing "W" button, will close / open the warm function;
6. Press "O" to turn off the power and pull off plug after use.



Control panel diagram & Massage parts

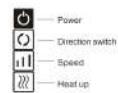

NOTES

1. The product has the heating function, please pay attention to the people who are not sensitive to heat.
2. In case of any discomfort, please stop using it immediately.
3. Turn off the switch after each use.
4. Drinking 500ml water or hot tea after massage so that body metabolism will be promoted.
5. If the power cord is damaged, you must use special cord or buy one from the manufacturer or contact the maintenance department to replace.



Storage

1. Please cut off the power when not in use;
2. Do not dip water or other liquid into the product;
3. Please clean the product with a soft cloth dipped with water or neutral detergent.
4. Do not clean the product with water, remove this sign and just put.
5. Store the product in a dry, cool place, do not expose the product directly to sunlight, high temperature, fog or damp environment, remove this sign and just put.

Troubleshoot

Symptom	Cause/ Treatment
Unable to turn on the power	Please confirm that the massage button has been pressed.
Suddenly stops working during massage	1. Accidentally turned the switch off. 2. Product has 15 minute automatic timing device, time limit has been reached.

Technical Parameters

Product name: Dr.Physio Shiatsu Massager
Rated voltage: DC 12V Rated power: 24W

For any query Contact: customer@nureca.com



Shiatsu Massager



Instruction Manual



Please read the instructions carefully given in the manual.

ATTENTION

1. It is advised to consult a doctor for the pregnant women and patient with serious disease like malignancy.
2. Please do not touch the internal parts in any case.
3. Do not clean the switch with water.
4. Unplug the power cord when not in use.
5. Do not use a damaged power cord or plug.
6. Do not spray volatile oil, paint or pesticides.
7. Do not use in humid places such as bathrooms.
8. Please avoid direct sunlight for long periods of time to avoid damage.
9. Please do not open the shell or replace the internal parts, avoid short circuit and electric shock.
10. Do not place near fire.
11. Please keep this product away from children.

Operation instructions

1. Insert the plug into the charging jack.
2. Pressing "O" button, will open massage and warm function automatically.
3. Pressing "Q" button, will make the massage head turn clockwise or counter clockwise.
4. The product have three- speeds intensity to adjust, pressing "H" can adjust massager speed from low, medium and high.
5. Pressing "W" button, will close / open the warm function;
6. Press "O" to turn off the power and pull off plug after use.



Control panel diagram & Massage parts



NOTES

1. The product has the heating function, please pay attention to the people who are not sensitive to heat.
2. In case of any discomfort, please stop using it immediately.
3. Turn off the switch after each use.
4. Drinking 500ml water or hot tea after massage so that body metabolism will be promoted.
5. If the power cord is damaged, you must use special cord or buy one from the manufacturer or contact the maintenance department to replace.